



AESTHETICS

Precision | Innovation | Care

PRODUCT PORTFOLIO

About us

Welcome to the dynamic partnership between **Pulse** where healthcare , where healthcare excellence meets cutting-edge aesthetics technology. Together, we are dedicated to redefining the boundaries of healthcare and beauty by blending medical expertise with innovative LASER solutions.

Our Mission



Quality and Innovation



Customer-Centric Approach



Wide Range of Products



Reliability and Trust



Advanced Aesthetic LASER Machines



Index

1. Wavelength Pro X Series

2. Pi-Code Series

3. Code Factor

4. Nodd

5. Hydra Frac

6. Cool Shape

7. Em Code

8. Bbl Super IPL+DPL

9. Visage

10. Excimer

11. Em-Ma

12. Skin Analyzer

13. Facecode

Wavelength Pro X

Product Description

A breakthrough hair removal platform that combines ultimate functionality with unparalleled comfort. Thanks to its adequate 12x16mm² spot size and advanced cooling system, treatments are now much faster and pain-free offering a significantly improved solution for patients and practitioners alike. Leverages the benefits of 4 combined wavelengths, with an improved patient experience and a business-oriented approach, creating a unique and result-driven new solution in the world of professional hair removal.



4 Wave technology
755nm + 808nm + 940nm + 1064nm



Humanized handle design
Spot size = 12 x 16mm²



Continuous cooling system



Safe & effective



Suitable for all skin types (I - VI).



Treatments




Hair
Removal



Skin
Rejuvenation

Theory

The Wavelength Pro X LASER system can emit LASER light of (755nm+808nm+940nm+1064nm) wavelength. It can be absorbed by the melanin located in the hair follicle so as to destroy the hair follicle. The hair shaft can also absorb the energy from LASER light, so as to be destroyed and gasified. Then the unwanted hair can be removed completely and permanently.

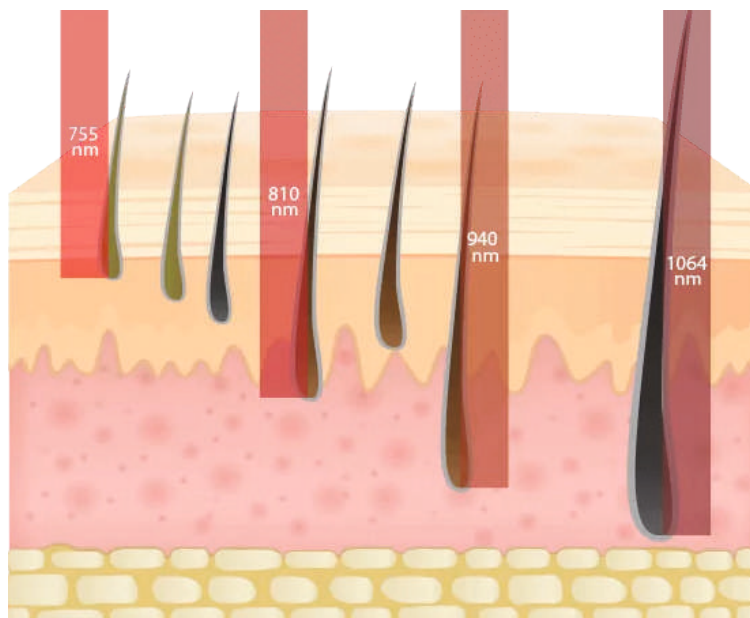

4 Wave Technology

755nm

Light and
Fine Hair Type

808nm

Effective & Safe
Wavelength

940nm

Deep & Coarse
Hair Type

1064nm

Deep & Fine
Hair Type



Technical Specification

TYPE OF LIGHT SOURCE	Diode LASER
WAVELENGTH	755nm+808nm+940nm+1064nm
PULSE WIDTH	5-400ms
FREQUENCY	1-10Hz
SPOT SIZE	12mm x 16mm
ENERGY DENSITY	Max 2-100J/Cm ² (Adjustable)
COOLING SYSTEM	Closed cycle water cooling system
OUTPUT POWER	1600W/2400W/600W (optional)
VOLTAGE	230VAC, 50-60Hz
NET	57 Kg/81 Kg
DISPLAY	15" Color touch LCD display

Pi-Code

Product Description

Pi-Code presents a new era in skin rejuvenation and tattoo removal, offering optimal results for wrinkle reduction, removal of colourful tattoos, and beginning pigmented lesions with just a few treatments.

Pi-Code mode of operation is based on delivering ultra-short picosecond pulses of energy into the skin.



A versatile platform enabling multiple patient-tailored treatment options.



Shortest pulse duration, enabling high peak power.




Unique double pulse mode



Supreme feature set



Suitable for all skin types.

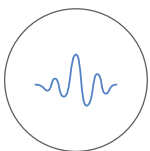


Double mode long pulse



Theory

Utilizing the explosive effect of the Nd: YAG LASER, the LASER permeate the epidermis into the dermis which includes amount of pigment mass. Since the LASER pulses in nanosecond but with super high energy, the shot pigment mass swell quickly and breaks into small pieces, which will be eliminated through metabolic system.

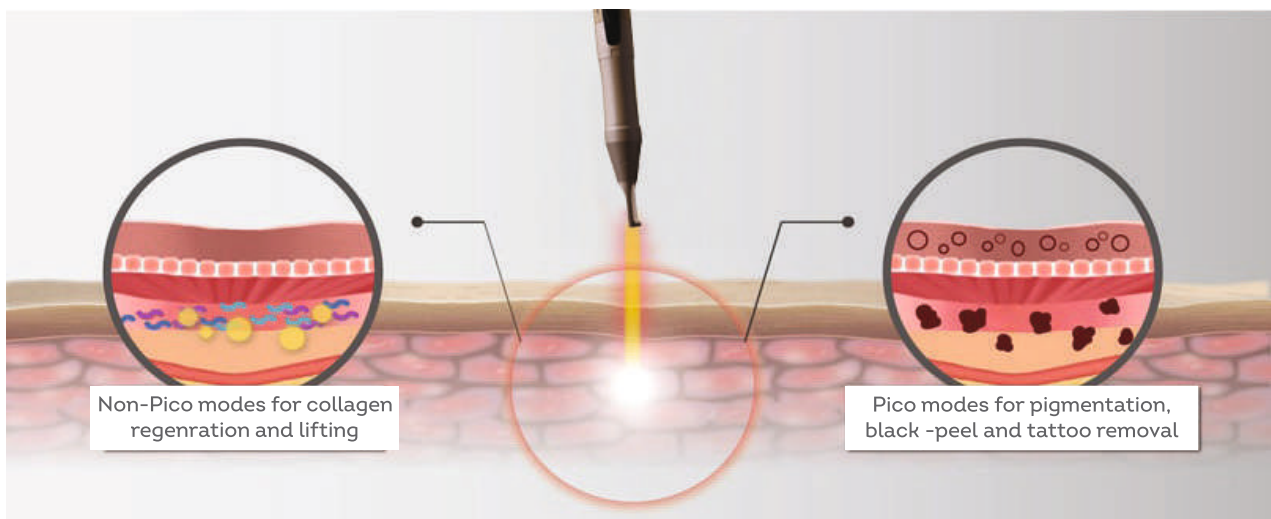


Treatments



Skin
Rejuvenation



Tattoo
Removal

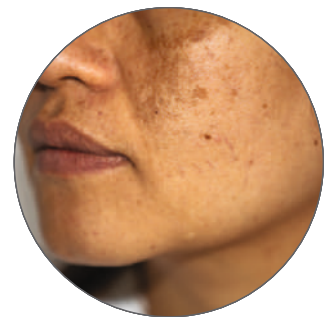

Pigmentation
Lesions

Technical Specification

SYSTEM TYPE	Picosecond	Long pulse
WAVELENGTH	1064nm 532nm	1064nm 532nm
PULSE WIDTH	2000ps	SLP : 300µs MLP :1-5000us
MAXIMUM REP RATE	Up to 10 Hz	Short: Up to 10 Hz Long: 1hz
SPOT SIZE	2-10mm (Single pulse: 1064nm: 100-800mj;532nm:100-400mj)	2-10mm
MAXIMUM ENERGY	Double pulse :1064nm:100-1600mj,532nm:100-800mj	SLP: 100-3000mj MLP: 1-30J
LIGHT GUIDE BEAM	650nm Red LASER	650nm red LASER
POWER REQUIREMENT	220V,10A,50/60Hz;	
MACHINE DIMENSIONS	680mm×400mm×820mm	
NET WEIGHT	80Kg	

Code Factor

Product Description


Code Factor CO₂ Fractional LASER, a state-of-the-art, minimally-invasive outpatient solution for treating stretch marks, scars, wrinkles, and skin dermabrasion. This cutting-edge device incorporates CO₂ technology to offer fast, painless, and safe procedures with highly effective results



360° Degree vaginal treatment option



40 Watt USA coherent RF excited laser source

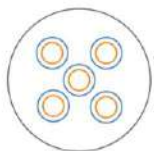

5 Modes to achieve optimum results



No consumables



Korean high quality 7 joint articulated arm



Treatments



Skin
Rejuvenation




Skin
Resurfacing



Post-Delivery
Rehabilitation



Scar
Remodeling

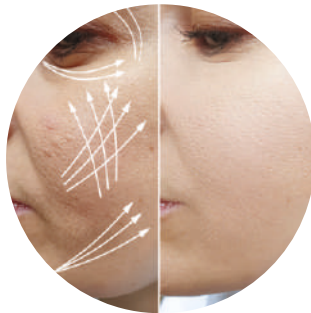

Skin
Tightening



Post-Menopause
-GSM



Vaginal
Tightening




Stress Urinary
Incontinence (SUI)



Vaginal Dryness &
Recurrent Infections

Fractional Mode

Seven Different Shapes

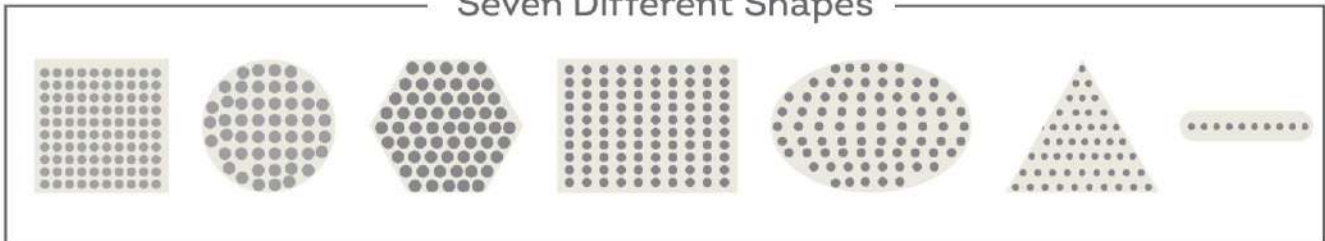

*Fractional CO₂ LASER
for vaginal tightening*



*Fractional CO₂ LASER
for vulva treatment*



*CO₂
for surgery*



*Fractional CO₂ LASER
for facial treatment*


Technical Specification

LASER	DC Closed RF LASER generator
WAVELENGTH	10.6um
WORKING MODES	Continuous/Ultra-pulse/Non-sequential pulse
SIZE	530*570*1380mm
POWER	40W
PULSE WIDTH	0.1-10ms adjustable Step 0.1ms
VOLTAGE	110v-240 AC
INDICATOR LIGHT	635nm infrared
LASER ORIGIN	USA
SCANNING GRAPHICS	Square, rectangle, circle, diamond, triangle, ellipse.
LANGUAGE	English
LCD DISPLAY	15 Inch 16- bit color LCD Display
COOLING SYSTEM	Air


Nodd

Product Description


Nodd Fractional RF Microneedle is designed specially for delivery RF energy directly to epidermis or dermis layer in order to stimulate and re-model collagen effectively by either RF micro-needling mode(invasive) or RF matrix mode (non-invasive). It can be used for both epidermal application and dermal applications in one platform without change any operational environment.

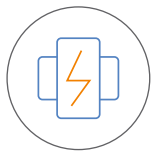


Stepping motor
type

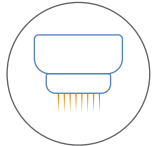

Golda plating
needling

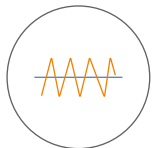

Mono-polar
bi-polar

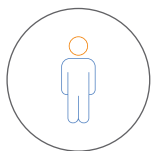

All body
treatment



Treatment



Fine
Lines



Deep
Wrinkles



Large
Pores



Acne
Scars



Stretch
Marks




Injury or
Surgery Scar




Loose &
Sagging Skin




Principles of Treatment

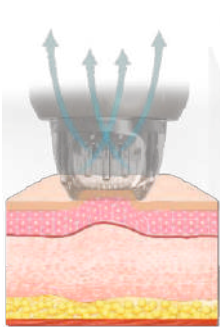


Auto Capture
Skin

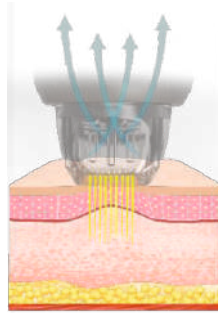

Creat
Channels

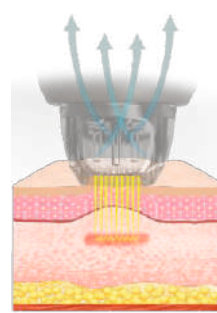

Release vacume
& Needle

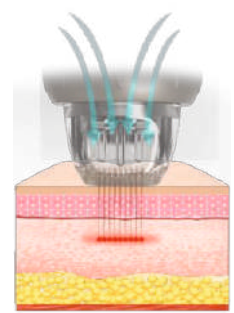

Finished


Micro Needle RF Applicator



M-18
Stepping
Motor Type



M-25
Needle
Depth Control

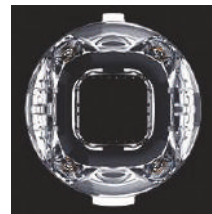

M-49
Gold
Plating



M-81
Safty Needle
System


Fractional RF Application



F-3*3
Eyes
Nose Area



F-5*5
Face
Area

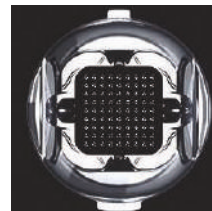

F-9*9
Body
Area

Technical Specification

POWER	Max 50 watts
HANDLE	Fractional RF handle micro needle RF handle
FREQUENCY	2MHz
ENERGY	1-10 Level
PULSE DURATION	10-900ms
NEEDLE DEPTHS	0.4-4mm Adjustable
PRESSURE LEVEL	1-4
WORKING MODE	Single/auto
MODE	Mono & bipolar
COOLING	1-3 Level
TIPS	Insulated & non-insulated
INVASIVE TIPS	M-18, M-25, M-49, M-81
NON INVASIVE TIPS	F-3*3, F.5*5, F-9*9

HYDRAFRAC

Product Description


Hydrafrac is an outstanding standalone treatment, as well as a must-have pre-treatment solution used to prepare skin for various skin rejuvenation treatments. It also offers amazing results as a pre-treatment for a diverse range of skin procedures.



Infusion of three technologies



Customised protocols

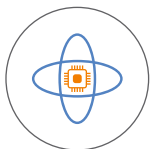


Universal applicability



Treatments



Skin
Hydration

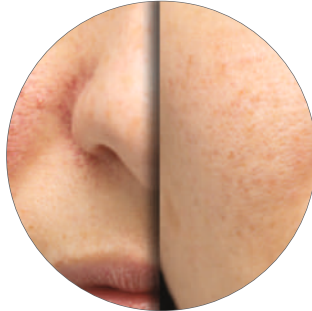

Skin
Rejuvenation



Deep
Cleansing



Blackheads
Removal




Wrinkle
Removal

Technology

Non-invasive microneedling apparatus



HydraFrac Aqua Peel

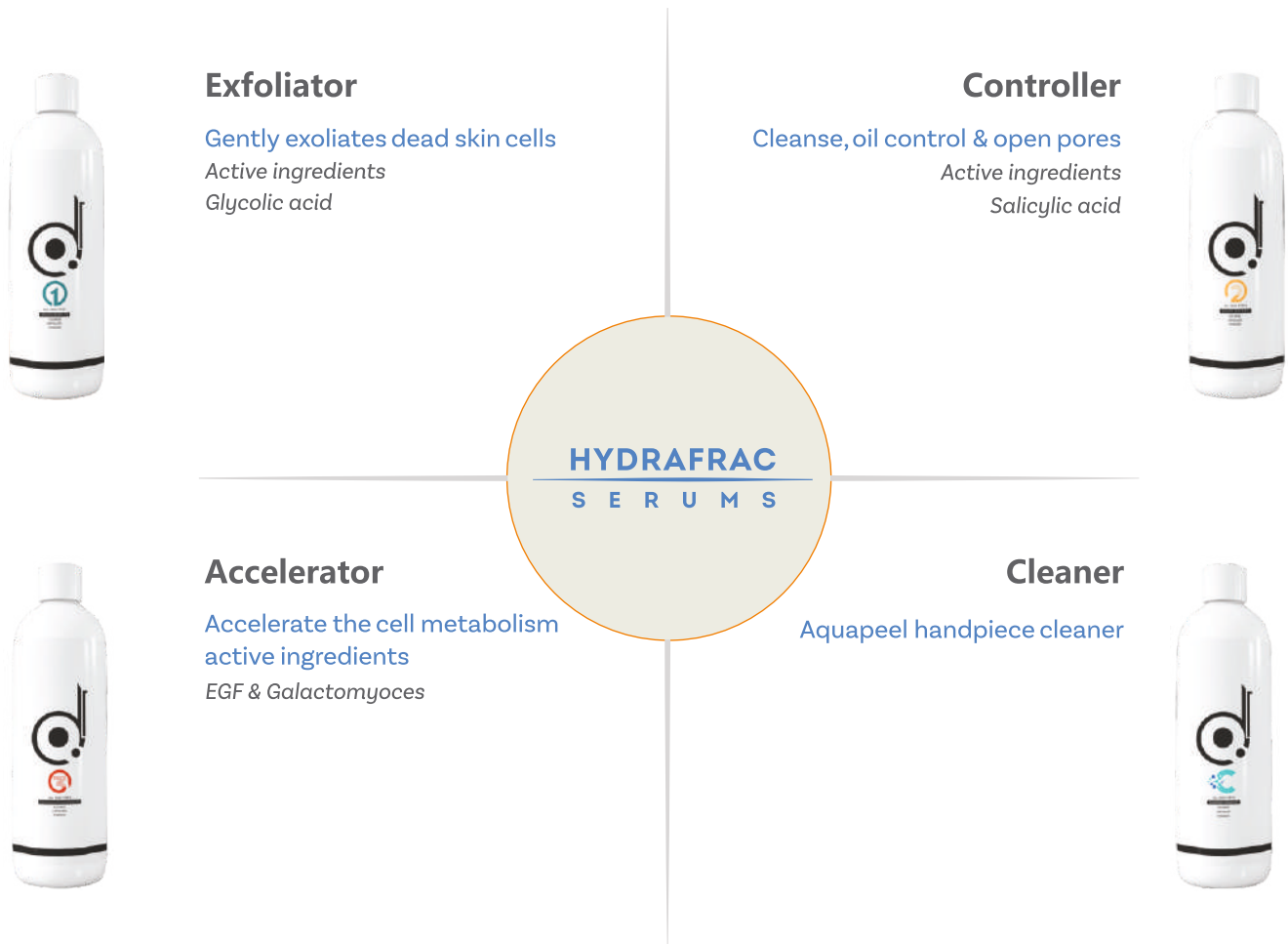

Hydrafrac Glow-o



Best for exfoliation and targeting pigmentation



Hydrafrac SONO

Best for product infusion and targeting dehydrated skin

Hydrafrac MRF

Best for skin-tightening



HTCT

Best for evening skin tone & invigorating the skin





Electro-Ion

Best skin lifting and uses electroporation

Technical Specification

SCREEN	10 inch TFT
ULTRASOUND	1Mhz/2ww/cm
INPUT VOLTAGE	110V/220V 50/60Hz
VACUUM RANGE	Max 680nm HG
ION LIFTING	500Hz (Digital Ion lifting)
RADIO FREQUENCY	1Mhz, multi-polar (6pol)
NOISE LEVEL	45DB

Cool Shape

Product Description


Cool shape uses state-of-the-art freezing technology to remove excess subcutaneous fat from a specific body part by dissolving fat cells using cool temperature without damaging the surrounding tissue. Once cold, fat cells will decompose naturally and gently be eliminated through the body's normal metabolic processes. designed with 360°-Degree cooling technology. The multiple cooling cups fit all curves and contours from as high as your chin down to the knees to deliver absolute, non-invasive body contouring your way.



360°
Cooling



Medical
silicone

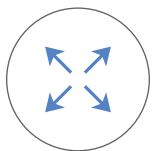

7 in 1 multi-size
applicator



Medical flexible LSR
silicone cups



Treatments



Body
Slimming



Double
Chin



Fat
Reducing

Theory

7 IN 1 MULTI-SIZE APPLICATOR

Ergonomic design, skin-friendly fit the body arc meet the whole body fat reduction needs from the face to the legs

Medical flexible LSR silicone cups

Whole shaped molding 360 cups dedicated and soft, comprehensive care for the skin.

360° Cooling


Better Cooling effect and less treatment time.

Liquid Silicone

Soft and tender touch like a baby skin.

Dual Temperature


Exclusive patent assures safety and efficacy.

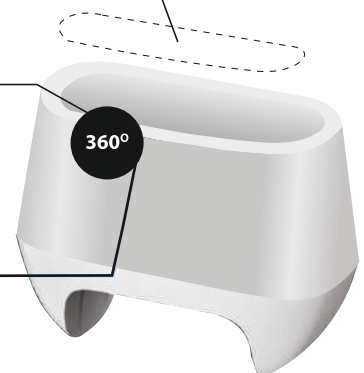


Flexible LSR Silicone Cups



Treatment Application

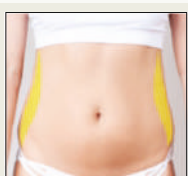

You may struggle with stubborn fat even after jogging and dieting.

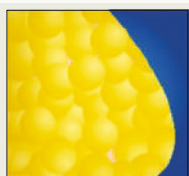


The coolslimming suction your skin and subcutaneous layer.




Cooling energy then freezes fat deposit to induce apoptosis.



Crystallized fat cells die off through the body's natural metabolism.

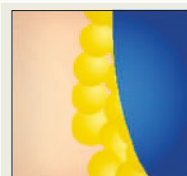

Remaining fat in the e target area became leaner and less dense.

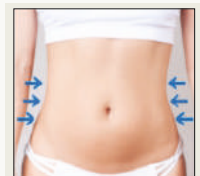

Combine with your healthy regimen to look slimmer and confident.

Technical Specification

BASIC OPERATION	
TYPE OF EQUIPMENT	Controlled cooling device for medical use
TECHNICAL CHARACTERISTICS	
MINIMUM TEMPERATURE	- 11°C
HAND-PIECES	Straight hp, curved hp, tight hp double hp, tiny hp, oval hp
POWER SUPPLY	100V-240V
MAXIMUM POWER	1000W
COOLING CAPACITY IN EACH HANDPIECE	90 W Programming $\pm 0^{\circ}\text{C}$ With a Temp of $\pm 25^{\circ}\text{C}$
REGULATION	
CLASS OF EQUIPMENT AS STANDARD	Class I, type BF
CLASSIFICATION ACCORDING TO DIRECTIVE	Class II a
THE DEVICE COMPLIES WITH THE REQUIREMENTS OF DIRECTIVES	Directive 2007/47 CE (93/42/CE) Of directives (medical device) directive 2002/95 CEE (RoHS) directive 2002/96 ECC (RAEE)
DIMENSIONS	540 x 440 x 550 mm
WEIGHT	39 kg
SCREEN	Touch screen 10,4"


Em Code

Product Description

Em Code is an effective way to improve your fitness level, body image and well-being in far less time. With the latest generation technology, high-tech non-water cooling system, high-quality insulated material, scientific magnetic pulses, and automatic system, one size is suitable for multiple areas. Em Code body workout activities over 90% of the muscle fibres at the same time and also deeper muscles are trained.



Optimised
applicator

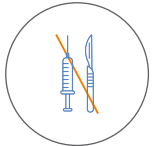

Non
invasive



Non
radiative



Treatments



Butt
Lifting




Muscle
Building



Fat
Reduction

User Friendly Screen

- High resolution interface
- Free sliding operation
- 45° Adjustable for best viewing
- Auto adjustment energy range when changing frequency
- Smart control of total tesla per session
- Safety low frequency waves without adverse effect

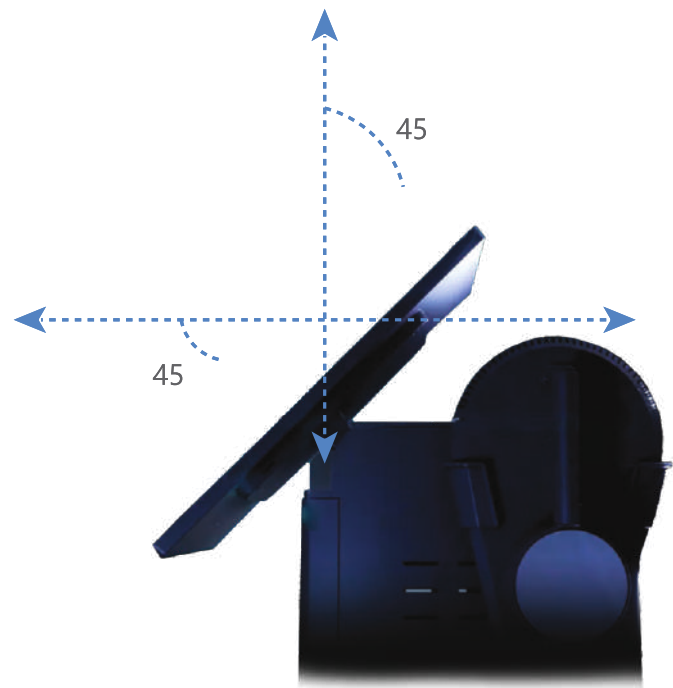

Optimized Applicator

- Latest generation technology
- High-tech non water cooling
- High quality insulated material
- One size suitable for multiple areas
- Unique up to 20,000 contraction per session more than any other equipment on the market



Technical Specification

PULSE FREQUENCY	1-150Hz
EM ENERGY	0-100%
RF ENERGY	0-100%
RF FREQUENCY	6MHz
PULSE TIME	285us
MAGNETIC FIELD STRENGTH	0.4 - 2.5T
TREATMENT TIME	Up to 60 min
EQUIPMENT WEIGHT	40Kg
USER INTERFACE	Touch screen
OPERATION TYPE	Continuous

Bbl Super IPL+DPL

Product Description

Intense Pulsed Light (abbreviated IPL) describes the use of intense pulses of non-coherent light distributed over a range of wavelengths from 400nm to 1200nm. The technology utilizes specific wavelength ranges to target various chromophores in the skin. This enables effective treatment of a broad spectrum of conditions, including vascular and pigmented lesions, skin photoaging treatments and hair removal.



Adjustable spot size design in handpiece



Hot and cold dual modes



Comfortable handpiece



Convenient to storage the filters




Universal handpiece with eleven filters



Treatment



Superficial
Pigment

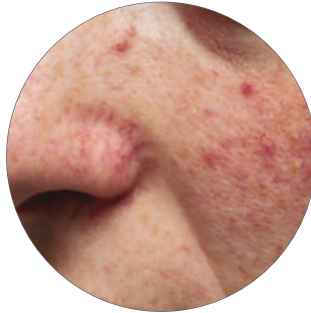

Superficial
Vessels



Pigmentation



Anti
Dullness



Whitening



Anti
phlogosis



Anti
Redness



Skin Fine
Lines



Hair
Removal

Precision Treatment Head (Optional)



38 X 10



10 X 10



φ8



φ6

Super Comfortable

Precise Operation, No Harm on Skin




Technical Specification

FILTERS	Standard filter	{ 420nm (420-950nm),515nm (515-950nm) 560nm (560-950nm),590nm (590-950nm) 615nm (615-950nm),640nm (640-950nm) 695nm (695-950nm)
	Optional filter	{ 500nm (500-600nm),550nm (550-650nm) ACNE (400-600 & 800-1200nm) VASCULAR (530-650 & 900-1200nm)
SPOT SIZE	14*38mm	
HANDPIECE	IPL & DPL handpiece calming & anti-sensitive handpiece (optional)	
FLUENCE	2-50J/cm2	
FREQUENCY	1-30Hz	
PULSE WIDTH	1-30ms	
PULSE WIDTH	1-3	
PULSE INTERVAL	Interval, 1-99ms	
SCREEN	10.1 Inch Industrial capacitive touch screen	
VOLTAGE	AC220V 50HZ	
POWER	2300W	
COOLING SYSTEM	Tec electronics cooling +water cooling+ air cooling	
WORKING TEMPERATURE	15-30 °C	
DIMENSION/G.W	Machine: 810x540x675mm, 28 KG Trolley: 645 515 1000mm, 20KG (Optional)	


Visage

Product Description

Visage targets & lifts frontalis muscles, corners of the mouth, zygomaticus major and minor muscles, matching the contours of every patient's face, activating & contracting these individual muscles. Meanwhile, special temperatures induce new collagen & elastin growth to resist aging. The combined effects restore and elevate the support of facial tissues, leading to an overall improvement in facial lifting.

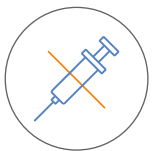


Needle-free
solution

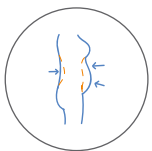

Simultaneous skin &
muscle enhancement



Natural
lift

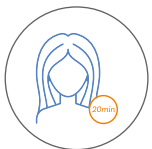

Full face in 20
minutes



Hand
free



Treatment



Wrinkles



Muscle
Tone



Facial
Lifting



Collagen

Simultaneous Skin & Muscle Enhancement



- *Injection-Free Solution*
- *Hands-Free*
- *Quick Results*
- *No Downtime*

- One Course Offers long-term results with the overall improvement of the structure of the face and skin quality
- No downtime required, nor the need for at-home recovery
- Gentle and pain-free non-invasive facial treatment requires no scalpels, no needles, no fillers
- It is safe for most people, including those with sensitive skin
- Equipped with unique technologies to target all facial layers, delivering non-invasive facial rejuvenation

Technical Specification

MAIN FREQUENCY	50-60Hz
INPUT VOLTAGE	100-220V
CHINA APPLICATOR WEIGHT	500g
CHEEK APPLICATOR WEIGHT	900g
MAX RF OUTPUT POWER	50W
RF FREQUENCY	1Mhz
DIMENSION	44*44*112cm
NET WEIGHT	30Kg

Excimer

Product Description

The Excimer 308 system is a handheld system for intense monochromatic UVB therapy, treating autoimmune diseases of the skin such as psoriasis and vitiligo. Using a 308nm wavelength which replicates the wavelength of the excimer laser, the excimer 308 system allows you to treat patients quickly and effectively with the benefits of a lighter-weight system that offers greater flexibility, lower investment costs, no consumables and greater ease of use.




High dose, highly effective and fast



No replacement consumables, low maintenance cost



User-friendly interface, convenient and easy to operate

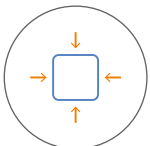

Compact, portable system



308nm Single wavelength



Treatments



Skin Lesions
and Vitiligo



Alopecia



Localised
Scleroderma

Theory

The Excimer 308 system uses a 308nm wavelength to effectively target skin lesions. The 16 cm² spot can be applied directly into localized plaques, or masks may be mounted to adapt the spot size and shape to smaller areas. A variable power density ranges from 50 to 6,000 mj/cm² and can be adapted to The needs of the individual patient. The system enables you to apply high, therapeutically effective doses of monochromatic UVB light while leaving healthy skin unexposed and protected.



Step-1

Vitiligo: localised,
few lesions.



Step-2

Psoriasis: localised,
palmoplantar, scalp
psoriasis.




Step-3

Alopecia areata




Step-4

Cutaneous T cell
lymphoma

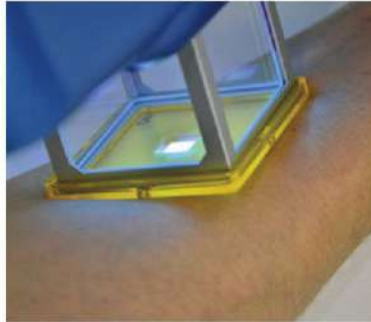

Step-5

Localised scleroderma

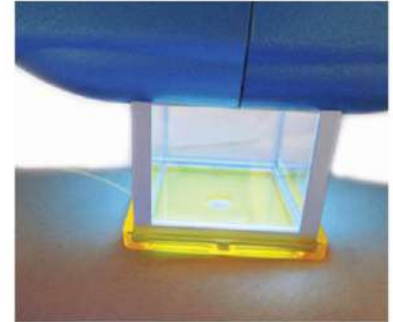

Step-6

Granuma annular


Technical Specification

LAMP TYPE	XeCl lamp
LIGHT WAVELENGTH	308nm
OPTICAL POWER DENSITY	308mW/cm ² +20%
FLUENCY	30-4200mJ/cm ²
SPOT SIZE	Max 16cm ²
TREATMENT TIME	1s-140s
WORK MODE	Continuous
DIMENSION(LXWXH)	27x24x26cm
POWER REQUIREMENT	100-240V.47-63Hz

Em-Ma

Product Description

Em-Ma utilizes High-Intensity Focused Electromagnetic technology (HIFEM) to cause deep pelvic floor muscles stimulation and restoration of the neuromuscular control. Key effectiveness is based on focused electromagnetic energy, in-depth penetration and stimulation of the entire pelvic floor area. A single Em-Ma session brings thousands of supramaximal pelvic floor muscle

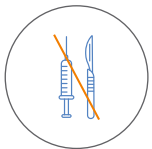

Non
invasive



Non
radiative



Patients remain
fully clothed



Thousands contractions
per session



Non
thermal



Operator
independent



Application


WOMEN



Pregnancy



Childbirth

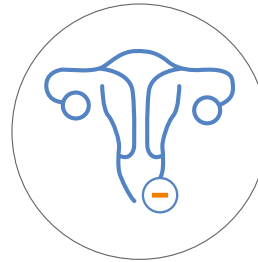

Hysterectomy



Menopause

MEN



Prostate
Problems



Bladder
Damage




Urinary
Stones



Sphincter
Problems


BOTH MEN & WOMEN

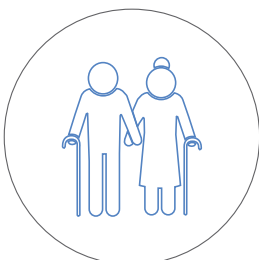


Aging



Neurological
Disorders

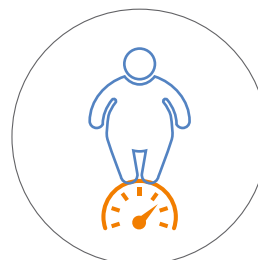


Obesity

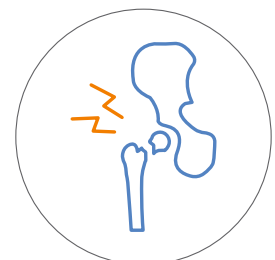


Pelvic
Trauma



Theory

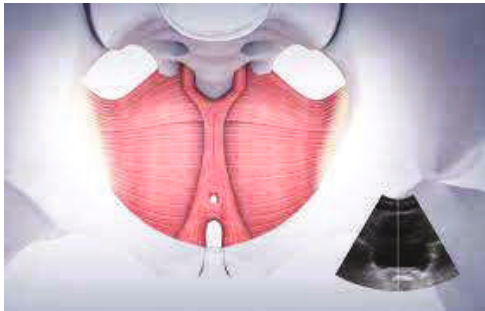


Myorram

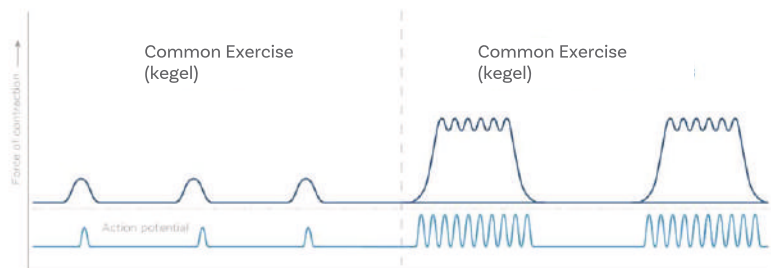

Em-Ma utilizes High-Intensity Focused Electromagnetic technology (HIFEM) to cause deep pelvic floor muscles stimulation and restoration of the neuromuscular control.


Backed by science

Em-Ma utilizes High-Intensity Focused Electromagnetic Technology (HIFEM) to cause deep pelvic floor muscles stimulation and restoration of the neuromuscular control.

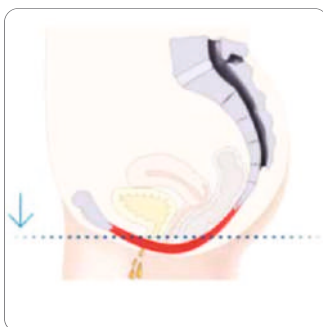

Before

- Pelvic floor muscles insufficiently support pelvic organs and affect bladder control.

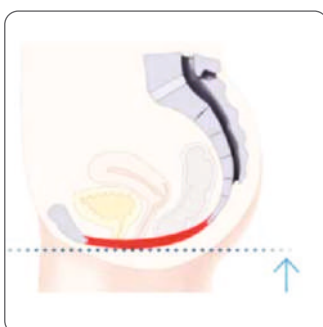


After

- Stimulation leads to regained control over pelvic floor muscles and bladder.

Skin Analyzer

Product Description

A skin analyzer is a device or tool designed to assess and analyze various aspects of the skin's condition. These devices use advanced technology to provide detailed information about the skin, helping individuals and skincare professionals make informed decisions about skincare routines and treatments.

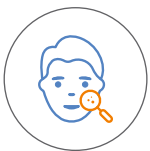

Skin digital analysis and quantification




Intelligent diagnosis report



8 Spectral imaging technology



36 Million pixel HD cameras

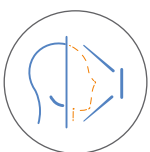

Automatic facial recognition



Interface Operation

8 SPECTRAL IMAGING TECHNOLOGY

RGB white light parallel polarized light cross polarized light blue light UV light Wushi light brown light compound light.

8 spectral imaging technology extracts skin images at different levels of the face for quantitative analysis.




Wu Spectrum



UV Spectrum



Cross Polarizedlight/Composite Spectrum

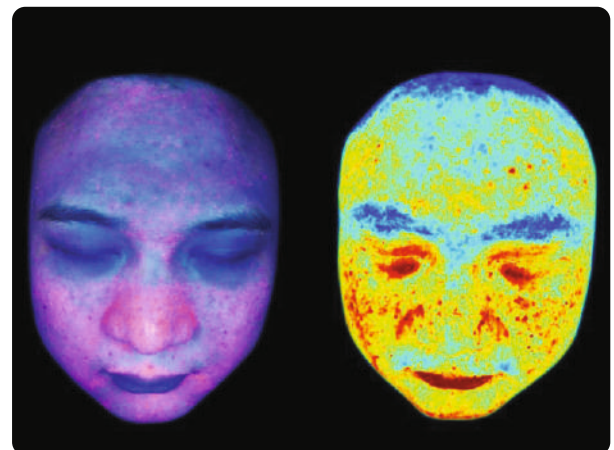

Blu-ray/Thermal Diagram

Technical Specification

Basic Operation

INSTRUMENT PARAMETERS	
OPERATING CONDITIONS: INDOOR U	System operating ambient temperature: 10-40
TRANSPORTATION AND STORAGE CONDITIONS	
TEMPERATURE RANGE: - 20° ~ 60°C	Humidity range: ≤ 80
TECHNICAL PARAMETERS	
POWER FREQUENCY: 50-60H	Power input voltage: global 100-240V



Facecode

Product Description

Facecode a non-invasive treatment, goes deeper to stimulate production of your body's own collagen and elastin, for natural-looking results without surgery or downtime.

Hand Piece

Thinner & curved line.
very comfortable to handle
double switches.

Transducers


Super thin

Display

LCD with touch screen

Main Body

Smiple & smart design

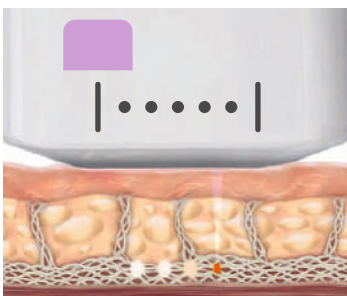

HIFU treatment Mechanism

HIFU is delivered into SMAS(Superficial Muscular Aponeurosis System), Fascia, deep dermis and substaneous fibrous tissue and generates the temperature increased, when the temprature grows more than 47°C (the Threshold of protein denaturation), the coagulation zone is generated. the zone is filled out by the tissue/muscle surrounded and in this moment, the skin muscle gets elastic and the wrinkle is tightened.

Application

Suitable for what kind of Conditions ?

Treat sagging skin, aging phenomenon for any areas such as forehead, eyebrow, cheekbones, apple muscles, cheeks, neck and body



Skin Tightening
and Rejuvenation



Non-Invasive
Lifting for Face



Body Contouring
and Tightening

Mutifunctional 2D handpiece



For Face Treatment



DS10-1.5
(Optional)



DS7-3.0
(Standard)



DS4-4.5
(Standard)

For Body Slimming




DS3-7.0
(Optional)



DS2-13.0
(Optional)

8D Handpiece

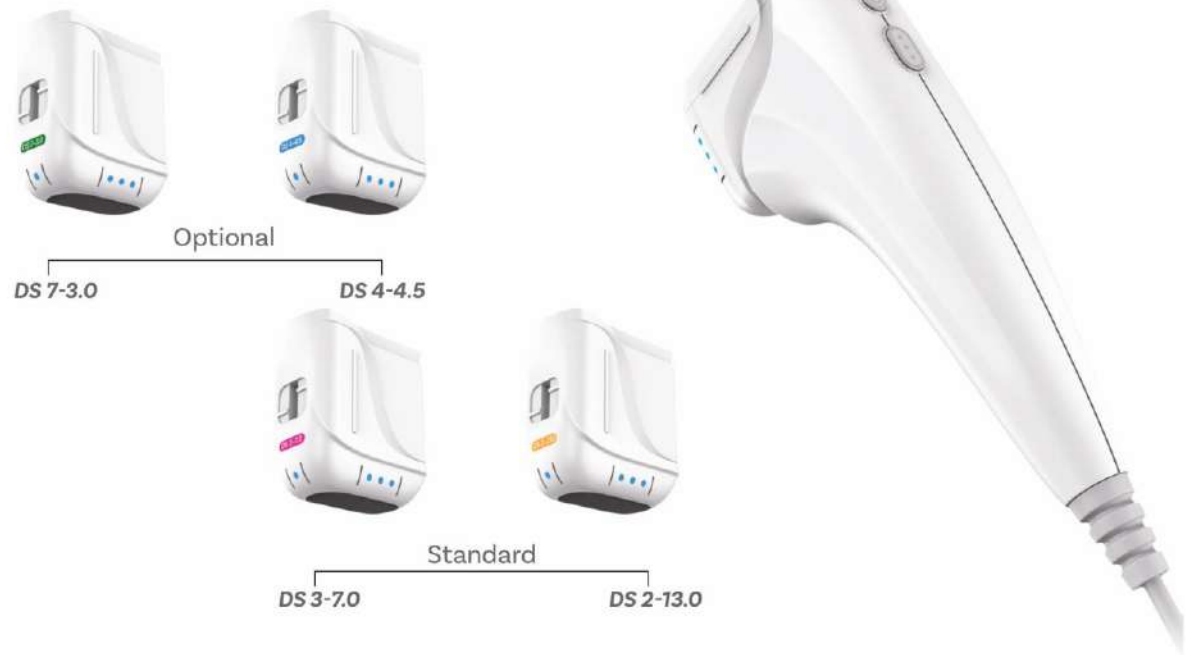

Technical Specification

POWER	Max 50 watts
HANDLE	Fractional RF handle micro needle RF handle
FREQUENCY	2Mhz
ENERGY	1-10 Level
PULSE DURATION	10-900ms
NEEDLE DEPTHS	0.4-4mm Adjustable
PRESSURE LEVEL	1-4
WORKING MODE	Single/auto
MODE	Mono & bipolar
COOLING	1-3 Level
TIPS	Insulated & non-insulated
INVASIVE TIPS	M-18, M-25, M-49, M-81
NON INVASIVE TIPS	F-3*3, F-5*5, F-9*9

ZuvioTechnologies Private Limited

1st Floor ' Plot no 646 , 27th Main Road 'HSR Layout ,
sec -1, Bengaluru 560102

P:+91 9899554722

E:info@pulseio.in

W:www.pulseio.in

